GB IE

BODY ANALYSER SCALE

Operation and Safety Notes

FI

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SE

VÅG SOM MÄTER KROPPSFETT

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
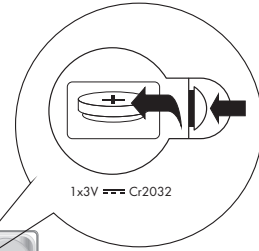
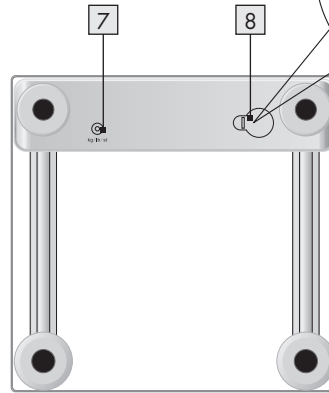
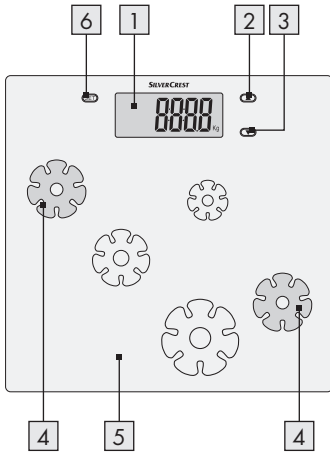
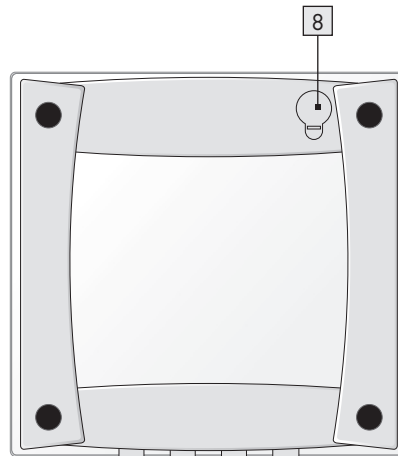
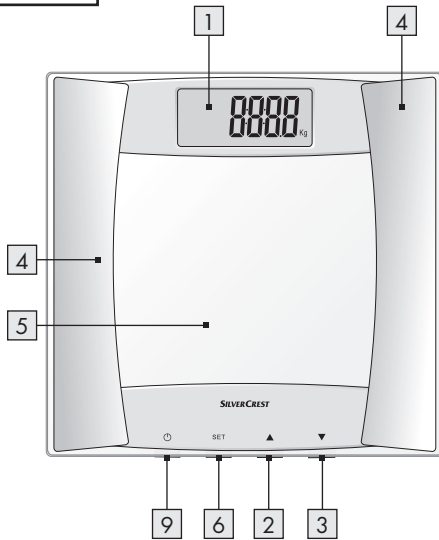
C

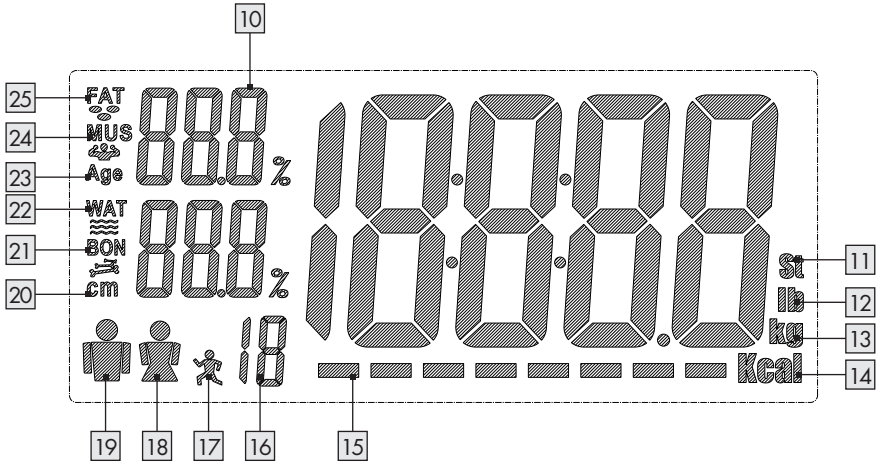Intended usage Page 6

Description of parts and features..... Page 6

Scope of Delivery Page 6

Technical Data Page 6

Safety notes Page 6

Before getting started

 Inserting / replacing the battery Page 8

Bringing into use

 Enter your personal data Page 8

Use

 Measuring the body weight Page 8

 Measuring the body weight, calorie requirement, proportion of body fat, water,
 bones and muscles..... Page 9

 Calibration of scales Page 9

 Meaning of the interpreted values..... Page 9

Trouble-shooting Page 9

Cleaning Page 10

Disposal Page 10

Body Analyser Scale

● Intended usage

These scales are intended for measuring the weight of people. The scales have a weighing range of up to 180 kg. Furthermore, the scales are designed to measure the body fat, water content, muscle and bone mass. Moreover, it is possible to determine the daily calorie requirement by entering the sex, age, height and weight. The device is not intended for commercial use.

● Description of parts and features

Scales

- 1 LC display
- 2 UP button
- 3 DOWN button
- 4 Electrode
- 5 Platform
- 6 SET button („settings“)
- 7 kg/lb/st button (weighing unit)
- 8 Battery compartment cover
- 9 ON/OFF button

Pictogram in the LC display

- 10 Percentage (of body fat, muscle, water, bone)
- 11 Unit of weight in „stones:pounds“
- 12 Unit of weight in „pounds“
- 13 Unit of weight in kilogram
- 14 Recommended daily calorie requirement
- 15 Body fat index
- 16 User ID (1-12)
- 17 Athlete mode
- 18 Female user
- 19 Male user
- 20 Body height in cm
- 21 Bone mass
- 22 Water proportion
- 23 Age
- 24 Muscle mass
- 25 Body fat proportion

● Scope of Delivery

Immediately after unpacking, check that all the items are present and that the device is in perfect condition.

- 1 Body Analyser Scale
- 1 Battery
- 1 Operating instructions

● Technical Data

Battery:	1 x 3 V lithium battery (CR2032 button cell)
Platform:	model Z31332 A/ Z31332B: 31 x 30 cm, safety glass
Load:	max. 180 kg (396 lbs)
Division (weight):	100 g (0.2 lb)
Division calories:	1 kcal
Memories:	12
Sex:	male / female
Age range:	6–100 years
Height range:	80–220 cm
Operating temperature:	10–40 °C

Weighing range:

Body fat:	3–50 %
Water:	25–75 %
Muscles:	25–75 %
Bones:	1–15 %

Safety notes

This safety advice is an important part of these operating instructions and must be observed every time electrical devices are used. The advice and instructions will help you to avoid risk of loss of life, accidents and damage. Therefore you must read carefully through the following information concerning safety and proper use.

- **⚠ WARNING! RISK OF FATAL INJURY AND RISK OF ACCIDENTS FOR INFANTS AND**

CHILDREN! Never leave children unsupervised with packaging materials or the device. Otherwise there is a risk to life from suffocation.

Always keep children away from the device.

- Children or other persons who do not have the knowledge or experience required to use the device, or who have limited physical, sensory or mental capabilities, must not use the device without supervision or instruction by a person who is responsible for their safety. Children must never be allowed to play with the device. The product is not a toy and must not be used as such.

⚠ WARNING! DANGER TO LIFE! Do not use the product if you have a heart pacemaker. Do not use the product in the vicinity of persons who have a heart pacemaker or use other medical devices.

⚠ CAUTION! RISK OF INJURY! Do not load the scales beyond 180 kg. The safety glass could break and cause cut injuries

⚠ CAUTION! RISK OF INJURY! Do not step on to the scales with wet feet and do not use the scales if the platform is moist. Otherwise you could slip.

- Protect your scales from impact, moisture, dust, chemicals and sudden changes in temperature and keep them away from sources of heat (ovens, heating radiators).
- Check the device for damage before bringing it into use. If the device is damaged do not use it.
- Remove the battery if the device breaks down or is defective. Otherwise there is an increased risk of leaking.
- Never open the housing of the device. This could result in personal injury or damage to the device. Always have any repairs carried out at the service centre or by an electrical equipment repair specialist.
- The product is not suitable for use by pregnant women.

Warning! Dangers arising from the use of batteries

⚠ WARNING! DANGER TO LIFE! Batteries are not intended for children. Do not leave batteries lying around. They may otherwise be swallowed by children or pets. If accidentally swallowed seek immediate medical attention.

CAUTION! EXPLOSION HAZARD! Never recharge single-use batteries. Do not short-circuit or open them. The battery may overheat, ignite or explode. Never throw batteries into fire or water. The batteries could explode.

Avoid extreme environmental conditions and temperatures as they may affect batteries e.g. above radiators.

Otherwise there is an increased risk of leaking. Exhausted or damaged batteries can cause caustic burns upon contact with the skin. To avoid this, always wear suitable protective gloves in these circumstances!

- Avoid contact with the skin, eyes and mucous membranes if fluid has leaked from the battery. In the event of contact with battery acid, thoroughly flush the affected area with clean water and seek immediate medical attention.
- Remove light dirt from the battery and product contacts with a dry, lint-free cloth before inserting the batteries.
- Remove the battery from the device if it is not going to be used for a prolonged period.
- When inserting the battery, ensure the correct polarity! This is shown in the battery compartment.
- Remove a spent battery immediately from the device. There is an increased risk of leaks!
- Batteries must not be disposed of in the normal domestic waste!
- Every consumer is legally obliged to dispose of batteries in the proper manner!
- Keep batteries away from children. Do not throw batteries into a fire. Never short-circuit them or take them apart.
- If the above instructions are not complied with, the battery may discharge itself beyond its end voltage. This will then pose a risk of leaking. If

the battery has leaked inside your device, you should remove it immediately in order to prevent damage to the device!

● Before getting started

Note: Hold the UP button [2] or the DOWN button [3] pressed down. This will speed up the process of setting the values.

● Inserting / replacing the battery

- The product is supplied with a battery already in place. Before the first use, pull the battery insulation strip out from between the battery and the battery compartment [8].

If the scales do not appear to work, check that the battery has been inserted correctly and reinsert it if necessary.

Note: The display indicates when it is time to replace the battery. When you step on to the scales and the battery is too weak, the LC display [1] shows „Lo“ and the scales switch off automatically. If you see this message, replace the exhausted battery with a new battery of the same type. When you insert the battery, make sure that the contacts do not get bent.

● Bringing into use

● Enter your personal data

1. For model A: Press the SET button [6] to switch on the scales. Press the SET button again to set your personal data.

For model B: Press the ON/OFF button [9] to switch on the scales. Then press the SET button [6] to set your personal data.

Note: The User ID [16] will flash in the LC display [1]. The number stands for the data profile of one person from the user group.

2. Press the UP button [2] or the DOWN button [3] to select the data profile (1-12).

3. Press the SET button to confirm the selected data profile as well as to set the option of the sex [18]/[19] and the athlete mode [17].

4. Press the UP button [2] or the DOWN button [3] to select your sex and the athlete mode.

Note: The athlete mode is a setting for sporty, highly trained or muscular persons.

5. Press the SET button [6] to confirm your entry.
6. Now make the further settings. Repeat steps 4 and 5 to set the age [23], height [20] and weighting unit [11]/[12]/[13].

After you have made the settings, the display will return automatically to the selection of the User ID. You can either make settings for further users or measure your body weight.

Note: The scales switch off automatically after about 10 seconds if they are not used.

● Use

● Measuring the body weight

Tip: Weigh yourself at the same time of day (preferably in the morning) after visiting the toilet, with an empty stomach and without clothing to ensure that your results are compatible.

- Place the scales down on a solid, flat horizontal surface.
- Step onto the platform [5] and do not move. Moving during the weighing process can lead to discrepancies in the results.

The scales will start to measure immediately.

After blinking twice, the weight display will show the result.

Note: The scales are overloaded if the display shows “Err”. Observe the maximum capacity of 180 kg / 396 lbs / 28 st.

- Step off the platform. The scales switch off automatically after about 20 seconds. Until then the weight measurement remains visible.

Note: The standard unit of weight is kilogram (kg). The factory default setting for the scales as supplied is “kg”. For model A, press the kg/lb/st button [7]

to choose between the units kg/lb/st:lb. For model B, press the SET button **6**.

● Measuring the body weight, calorie requirement, proportion of body fat, water, bones and muscles

- For model A: If the scales are switched off, press the SET button to switch on the scales.
- For model B: If the scales are switched off, press the ON/OFF button **9** to switch on the scales.
- Press the UP button **2** or the DOWN button **3** to proceed to your user profile.

Note: As soon as you reach your user profile, all your personal data will be shown in the LC display **1**.

Note: Press the SET button to adapt your personal data.

- Step barefoot onto the platform **5** and do not move. Ensure that you are standing on the two electrodes **4**.

The weight measurement is beginning; remain standing on the two electrodes. The measurement lasts a few seconds (display indication ---).

The LC display then shows the measured values three times alternating between two groups.

Group 1: Percentage of body fat **25** and water **22**, weight (kg/lb/st:lb) **11**, **12**, **13** and body fat index **15**

Group 2: Percentage of muscles **24**, bones **21** and calorie requirement **14**

Note: The display will show Err2, if you do not stand on the two electrodes during weighing, or if the total weight is above the weighing range.

The scales switch off automatically after these values have been displayed.

● Calibration of scales

Note: The scales must be re-calibrated whenever they are moved or after every battery replacement. To do this, proceed as follows:

- Step on to the scales. The scales will show you an incorrect measuring result.
- Step off the platform **5** once again. The scales will calibrate automatically. Wait until the scales switch off automatically.
- Step onto the scales once again. The scales will show you the correct measuring result.

● Meaning of the interpreted values

LC display	— —	— — — —	— — — — —	— — — — — — — —
Meaning	thin	normal	overweight	extremely overweight

● Trouble-shooting

The device contains delicate electronic components. This means that if it is placed near an object that transmits radio signals, it could cause interference. If the display indicates a problem, move such objects away from it.

Electrostatic discharges can lead to malfunctions. In cases of the appliance failing to work, remove the batteries for a short while and then replace them.

● Cleaning

- Do not use chemical cleaners or scouring agents for cleaning the device.
- Never immerse the device in water or other liquids. Ignoring this advice may damage the device.
- Clean the device using a slightly moist cloth. For more stubborn dirt you can also apply a weak soap solution to the cloth.

EMC CE

● Disposal

The packaging is made from environmentally-friendly materials, which may be disposed through your local recycling facilities.

Contact your local refuse disposal authority for more details on how to dispose of your worn-out product.

To help protect the environment, please dispose of the product properly when it has reached the end of its useful life and not in the household waste. Information on collection points and their opening hours can be obtained from your local authority.

Defective or used batteries must be recycled in accordance with Directive 2006/66/EC. Return batteries and /or the device to a local recycling facility.

Improper disposal of batteries can harm the environment!

Never dispose of batteries in your household waste. They may contain toxic heavy metals and are subject to hazardous waste regulations. The chemical symbols of the heavy metals are: Cd = cadmium, Hg = mercury, Pb = lead. Therefore dispose of used batteries through your municipal collection site.

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© by **ORFGEN Marketing**
OWIM GmbH & Co. KG